

Dominant 50



POWERFUL PERFORMANCE



High Vacuum Suction Pump

Dominant 50

High suction capacity and powerful performance

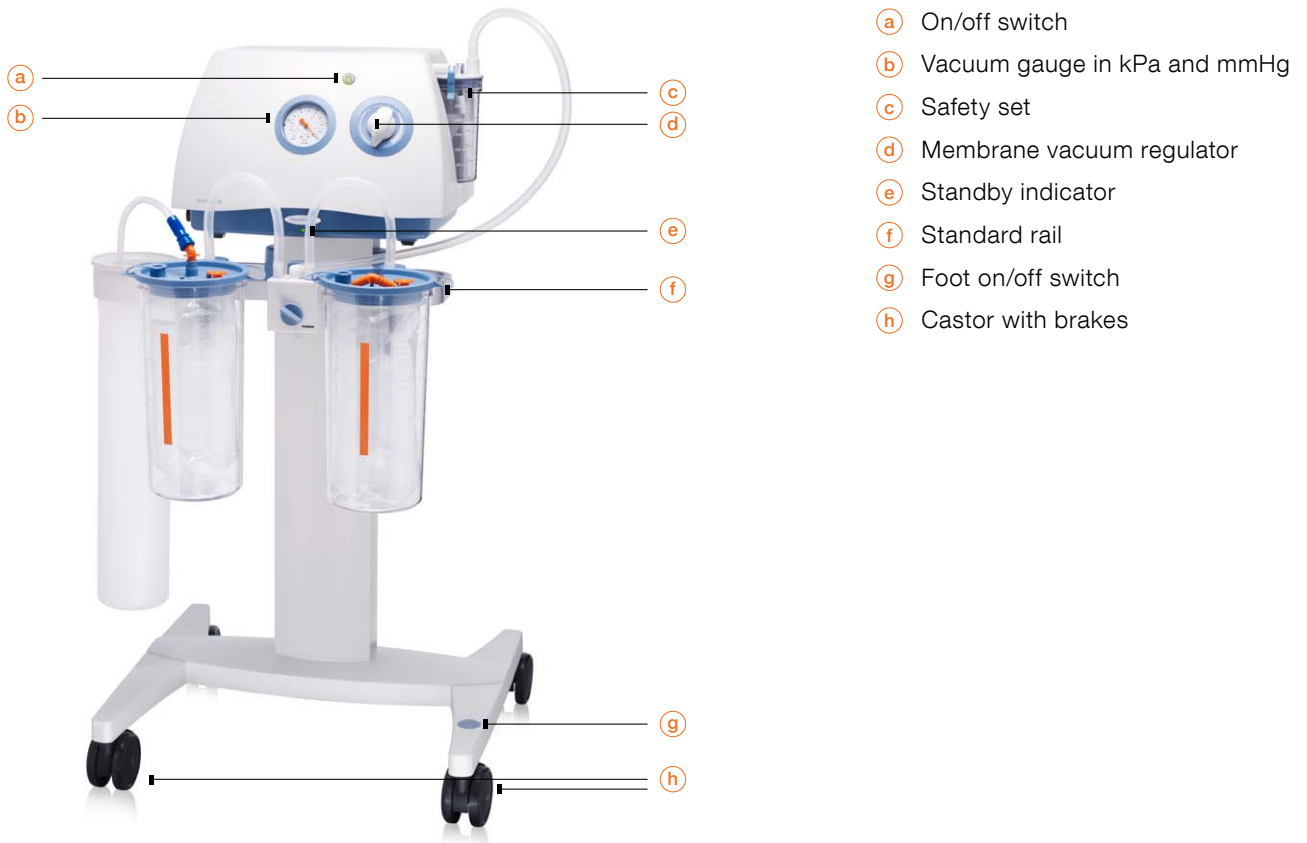
Fast vacuum build-up and strong performance

The **Dominant 50** is a powerful, high quality suction pump with the high suction capacity necessary for surgical suction and liposuction. In a range of settings, including the operating theatre and emergency care, the Dominant 50 can be adapted to suit to your specific and individual needs with a multitude of Medela accessories.

Medela's reliable piston/cylinder technology ensures a rapid vacuum build-up and delivers a high suction capacity. Thanks to the slow speed of rotation, Dominant 50 is vibration-free and amazingly unobtrusive, even when working at highest vacuum level. The membrane vacuum regulator allows for precise vacuum setting.

The **Dominant 50** is equipped with a foot on/off switch on the trolley, conveniently and hygienically leaving your hands free for surgery.

The **Dominant 50** has been awarded the prestigious Red Dot Award for outstanding design in the category Medicine and Rehabilitation.



Powerful and high suction capacity for a variety of surgical applications

The Dominant 50 reliably delivers efficient vacuum during surgery. Medela's piston/cylinder technology allows the vacuum build-up within seconds, which makes the Dominant 50 suitable for applications that require rapid suctioning of large quantities of fluids. The membrane vacuum regulator enables a precise vacuum adjustment.



Simple adjustment of vacuum level

Hygienic use and easy handling

The convincing design of the suction pump makes cleaning easy and hygienic. The safety set at the side of the pump prevents an overflow and reduces the risk of contamination. You can choose from a range of filters to protect the pump. A convenient foot on/off switch on the Dominant 50's trolley makes it easy and hygienic to operate during surgery.



Hygienic and easy to use

Flexible and adaptable to your individual and specific needs

The Dominant 50 is available as a rack and a mobile version and can be combined with a wide range of accessories including the Disposable or Reusable Collection System, foot vacuum regulator, clamp holders, tubings and a range of filters. The accessories can easily be attached to the rail of the pump and are designed to work perfectly together, offering the user a perfect complete system.



Hands free on/off switch

Fast Facts

- Medela offers healthcare professionals 50 years of expertise in medical vacuum technology.
- Our swiss development and manufacturing process yields mechanical precision to guarantee excellent performance.
- The fast vacuum build-up allows for efficient suctioning of larger amounts of fluids.
- The Dominant 50 is a high quality, long lasting suction pump, appreciated by medical technicians.
- The convincing design of the pump allows hygienic use and easy handling.
- We offer you a complete system, including Reusable and Disposable Collection Systems and accessories which are adaptable to your specific and individual needs.
- 5 year warranty.

Accessories

Foot vacuum regulator for hands free regulation of vacuum

can be attached to standard Medela suction jars and can be fixed to the standard rail when not in use.



On/off foot switch

for turning the Dominant 50 on/off. Equipped with a cable length of 3.5 m.



Various filters

to protect the pump from overflow, to protect the environment from bacteria or to neutralise odour.



Reusable Collection System

auto-clavable PSU jars and lids collect secretions easily and economically. Sizes: 1, 2, 3 and 5 litres.



Disposable Collection System

disposable suction liners and reusable PC jars offer a hygienic, simple handling and efficient fluid management.



Various tubings

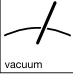
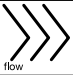

made from silicone or PVC.

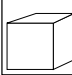



Sizes: 1.5 and 2.5 litres.

Use Medela accessories for even greater flexibility and ease of use.

Technical data

	high vacuum - 90 kPa - 675 mmHg
	high flow 50 L/min.
	17 kg 37.4 lbs

	H x W x D 965 x 510 x 480 mm 38 x 20 x 19 inches
	Power 100 - 120V 110W 230 - 240V 110W

	ISO 9001 ISO 13485 CE (93/42/EEC), Ila 5 years warranty	
		

Medical Vacuum Technology for Healthcare Professionals

Please contact us or your local Medela representative for details.

 Medela AG
Lättichstrasse 4b
6341 Baar / Switzerland
www.medela.com

Local contact:

Medela Australia Pty Ltd
Medical Technology
3 Arco Lane
Heatherton, Vic, 3202
Australia
Phone +61 (0) 3 9552 8600
Fax +61 (0) 3 9552 8699
contact@medela.com.au
www.medela.com.au

Medela Canada Inc.
4160 Sladeview Cres., Unit #8
Mississauga, Ontario
L5L 0A1
Canada
Toll Free 800.435.8316
Toll Free Fax 800.995.7867
info@medela.ca
www.medela.ca

Medela UK
Huntsman Drive
Northbank Industrial Park
Irlam, Manchester M44 5EG
UK
Phone +44 870 950 5994
Fax +44 870 389 2233
info@medela.co.uk
www.medela.co.uk

Medela Inc. USA
P.O. Box 660
1101 Corporate Drive
McHenry, IL 60050
USA
Phone +1 877 735 1626
Fax +1 815 363 2487
suction@medela.com
www.medelasuction.com